

Tracy, CA 95304
Phone: (866) 592-0488 • FAX: (209) 834-8832
topshelforthopedics.com



Stable Sling Universe

Instructions For Use

Please read these instructions completely and carefully before applying this brace. Correct application of the brace is important for proper function.








Stable Sling Universe

Instructions For Use

Please read these instructions completely and carefully before applying this brace. Correct application of the brace is important for proper function.

Application Instructions:

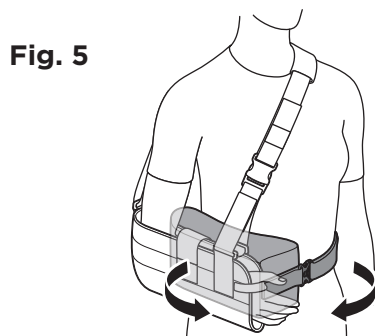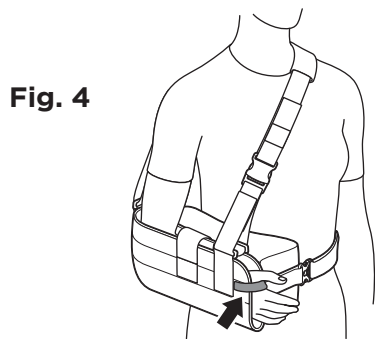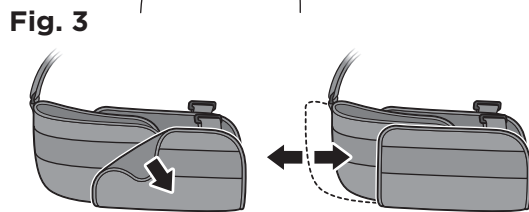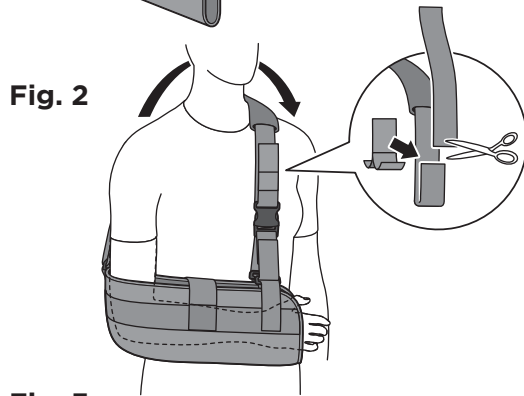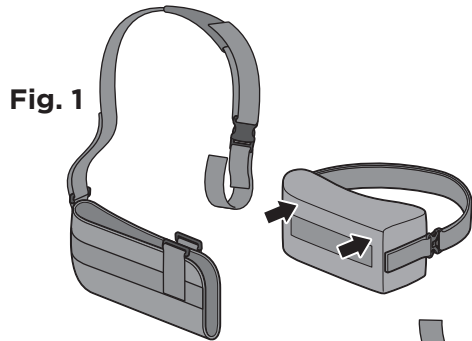

Indications: Rotator cuff repair, torn labral injuries, dislocation/subluxation of the shoulder, surgical repairs of the shoulder, shoulder strains, and capsular tightening.

Figure 1

Remove the sling from the pillow and detach the shoulder strap from the front of the sling.

Figure 2

Position the sling on the arm, making sure the elbow is all the way back in the sling. Secure the shoulder strap over the opposite shoulder to the front of the sling. Medical professional will adjust the length of the shoulder strap. Extra length may be trimmed by removing the Y-Tab Velcro piece, cutting off the extra, and then re-applying the Y-Tab Velcro. Minor adjustments may be made to the front of the strap below the buckle.

Figure 3

The sling should stop at the base of the little finger. If adjustment is required, the front piece may be moved forward or back to accommodate different forearm lengths.

Figure 4

If desired the small strap may be used in between the thumb and index finger to prevent the arm from sliding out of the sling.

Figure 5

Place the pillow at the waist of the affected side with the widest part of the pillow forward, secure the waist strap.

Use & Care: Hand wash in cold water using mild soap, rinse thoroughly.

Warning: This product is to be used under the supervision of a medical professional. This device is not intended for unsupervised public use. If you experience any pain, swelling, sensation changes, or any unusual reactions while using this product, consult your medical professional immediately.

Warranty: Top Shelf Orthopedics guarantees this product to be free of defects in material and construction for a period of six months from the date of purchase.

Caution: Federal Law (U.S.A) restricts this device for sale by or on the order of a licensed healthcare professional.